

## Data sheet


## Peltier-Cooled Incubator ILP 240 Smart



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

### Advantages of the SMART controller:

- 4,3", clear, full colour touch screen
- LAN, USB ports for data transfer
- multi-segment time and temperature programs
- visual and sound alarm
- internal memory for programs and data storage
- event registry
- user manual for direct download
- Quick change of program parameters
- Alarm Bar
- operating with gloves on

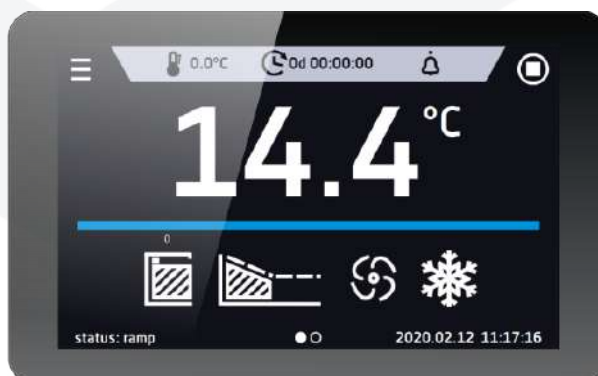

**Smart** - preview screen

## TECHNICAL DATA

air convection	forced
chamber capacity [l]	245
working capacity [l]	245
controller	microprocessor PID
display	4,3" full colour touch screen

## TEMPERATURE


temperature range [°C]	0...+70 (max 20°C below ambient temperature)
temperature resolution every ... [°C]	0,1
temperature fluctuation at 37°C [±/°C]*	0,1
temperature variation at 37°C [±/°C]*	0,3
temperature protection	class 2.0 to DIN 12880 / class 3.3 (option)

## CHAMBER

door type	double(5) / door with viewing window (option)
<b>interior</b>	
Smart	acid-proof stainless steel to DIN 1.4301
IG Smart	acid-proof stainless steel to DIN 1.4301
<b>housing</b>	
Smart	powder coated sheet
IG Smart	stainless steel linen finish

## overall dims [mm] /1/

width A	820
height B	1140
depth C	840
<b>internal dims [mm]</b>	
width D	600
height E	800
depth F	510

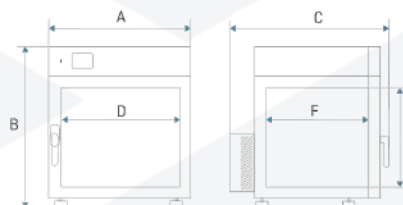

shelves (standard   max)	3   10
max shelf workload [kg] /2/	25
max unit workload [kg]	90
weight [kg]	140

## ELECTRICAL PARAMETERS

voltage**	230V 50-60Hz
nominal power [W]	800
warranty	24 months
manufacturer	POL-EKO

all the above technical data refer to standard units (without optional accessories)

\* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

$K = \pm (T \text{ average max.} - T \text{ average min.}) / 2$

\*\* - other power supplies on request


1 - depth doesn't include 50 mm of power cable, the width does not include the 20 mm of rubber plug

2 - on uniformly loaded surface

3 - reinforced shelf

4 - additional internal glass door

## OPTIONS AND ACCESSORIES

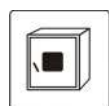

Order number: \*/A

door with viewing window



Order number: \*/P INOX

Stainless steel wire shelf INOX




Order number: \*/PP

Perforated shelf



Order number: \*/PW

Reinforced shelf

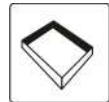

Order number: KUW GN\*/\*

Stainless steel cuvettes



Order number: QLK\*

Castors



Order number: \*/W

Reinforced version



Order number: LabDesk

LabDesk software



Order number: BRT\*/L or IQ/OQ/PQ

Calibration and IQ, OQ, PQ qualification



Order number: \*/3.3

Over temperature protection 3.3 class according to DIN 12880




Order number: BPP 12

Battery backup for display



Order number: PORT ALARM

Dry alarm contact

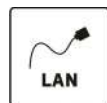

Order number: LANK

LAN cable



Order number: OCZ/20

Non-standard access port 20 mm




Order number: OCZ/30

Additional access port 30 mm



Order number: OCZ/60

Non-standard access port 60 mm

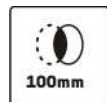

Order number: OCZ/100

Non-standard access port 100 mm



Order number: KD

Access control