

Data sheet


Ultra-Low Freezer ZLN-UT 200 VIP Smart



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

Advantages of the SMART controller:

- 4,3", clear, full colour touch screen
- LAN, USB ports for data transfer
- multi-segment time and temperature programs
- visual and sound alarm
- internal memory for programs and data storage
- event registry
- user manual for direct download
- Quick change of program parameters
- Alarm Bar
- operating with gloves on

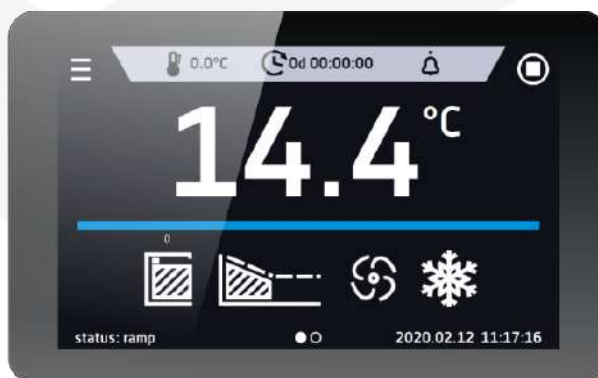

Smart - preview screen

TECHNICAL DATA

air convection	natural
chamber capacity [l]	259
number of boxes 133x133x50mm (option) [pcs]	192
controller	microprocessor PID
display	4,3" full colour touch screen
insulation type	vacuum - VIP panels

TEMPERATURE

temperature range [°C]	-86...-50
temperature resolution every ... [°C]	0,1
cooling time from +20°C to -80°C [min]	160
heating time in case of power failure -80°C to -60°C [min]	50
temperature fluctuation at -80°C [±/°C]*	1,5
temperature variation at -80°C [±/°C]*	4,0

CHAMBER


door type	solid
interior	
C Smart	stainless steel to DIN 1.4016
P Smart	acid-proof stainless steel to DIN 1.4301
housing	
C Smart	powder coated sheet
P Smart	powder coated sheet

overall dims [mm] /1/

width A	880
height B	1390
depth C	960

internal dims [mm]

width D	620
height E	770
depth F	580
height G	360

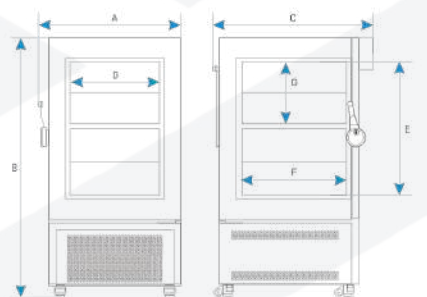

shelves (standard max)	2/2
max shelf workload [kg] /2/	10
number of internal chambers	2
max unit workload [kg]	65
weight [kg]	200

ELECTRICAL PARAMETERS

voltage**	230V 50-60Hz
nominal power [W]	2100
refrigerant	R290 / GWP=3 R170 / GWP=6
warranty	24 months
manufacturer	POL-EKO

all the above technical data refer to standard units (without optional accessories)

* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

$K = \pm (T \text{ average max.} - T \text{ average min.}) / 2$

** - other power supplies on request

1 - depth doesn't include 50 mm of power cable, the width does not include the 20 mm of rubber plug

OPTIONS AND ACCESSORIES




Order number LabDesk

LabDesk software



Order number: PORT ALARM

Dry alarm contact

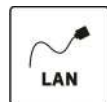

Order number: LANK

LAN cable



Order number: KD

Access control



Order number: ZLN-UT/CO2

CO₂ back up system



Order number: ZLN-UT/STP12 / 16

Boxes for samples



Order number: ZLN-UT/ST12 / 16

Rack with drawers